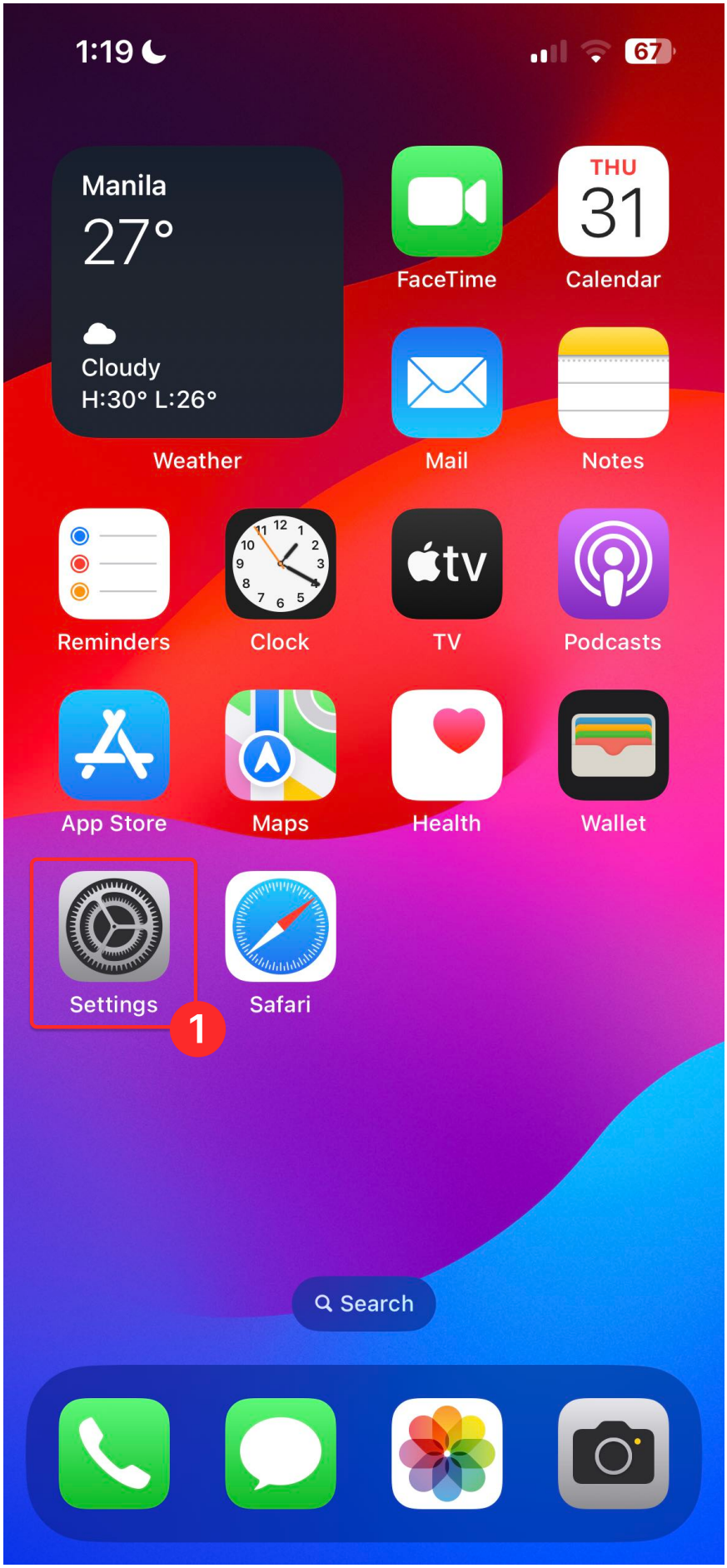
Stethoscope Mode in VSee Messenger (iOS)

Chi - 2024-10-30 - Telemedicine Equipment

On iOS devices, the stethoscope audio can be switched ON/OFF through the device's settings:

1. Click on the **Gear icon/Settings**.


1:19
67
Manila
27°
Cloudy
H:30° L:26°
Weather
FaceTime
THU
31
Calendar
Mail
Notes
Reminders
Clock
tv
TV
Podcasts
App Store
Maps
Health
Wallet
Settings
1
Safari
Search

1. Under Settings, scroll down and select **VSee Messenger.**



1. Scroll down. Under Advanced, toggle ON - **External Mic is Stethoscope**.

Note: Switch ON only when a stethoscope is connected to the device.



- Tags
- stethoscope
- stethoscope mode